



Oxford Cambridge and RSA

**Thursday 23 June 2022 – Morning**

**A Level Chemistry B (Salters)**

**H433/03** Practical skills in chemistry

Practical Insert

**Time allowed: 1 hour 30 minutes**



**INSTRUCTIONS**

- Do **not** send this Insert for marking. Keep it in the centre or recycle it.

**INFORMATION**

- This document has **4** pages.

### Determination of the percentage of copper present in brass.

Brass is an alloy of copper and zinc. The proportions of each metal depend on what the brass is being used for.


A group of students use colorimetry to determine the percentage of copper present in a sample of a brass.

The students' description of their method along with their results is shown below.

#### Constructing a calibration graph

1. A variety of standard copper(II) nitrate solutions (see results table) were made from an initial standard solution of concentration  $0.400 \text{ mol dm}^{-3}$ .
2. An orange filter was placed in the colorimeter.
3. A cuvette was rinsed and  $\frac{3}{4}$  filled with deionised water. The outside of the cuvette (see diagram below) was wiped dry with a paper tissue and placed in the colorimeter. The absorbance reading was adjusted until it read zero.
4. The deionised water was discarded and the cuvette  $\frac{3}{4}$  filled with the first solution of known concentration ( $0.025 \text{ mol dm}^{-3}$ ). The absorbance was measured and recorded in the results table.
5. The cuvette was again rinsed with deionised water and step 4 repeated for all the standard solutions. Their absorbance values were recorded in the results table.

A graph of absorbance against concentration can now be constructed.

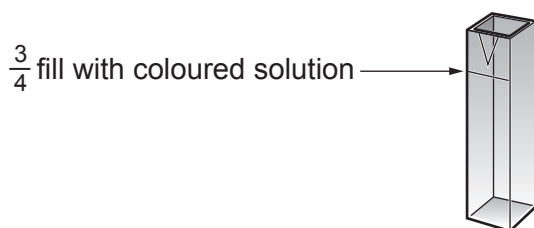

The cuvette that holds the solution in the colorimeter

### Preparation and measurement of the absorbance of a solution of the brass sample

1. Approximately 4 g of the brass sample were weighed accurately.
2. This mass of brass was dissolved in a minimum amount of concentrated nitric acid in a conical flask in a fume cupboard and the solution allowed to cool. The resultant solution was then carefully transferred to a 250 cm<sup>3</sup> volumetric flask and made up to the mark with deionised water.
3. Using a teat pipette, the cuvette used for constructing the calibration graph was rinsed with the brass solution,  $\frac{3}{4}$  filled, and its absorbance measured.  
This was compared with the standard solutions to determine the copper concentration.

### Results

Concentration of copper(II) nitrate solution / mol dm <sup>-3</sup>	Absorbance
0.025	0.05
0.050	0.12
0.100	0.25
0.200	0.40
0.300	0.72
0.400	1.03

Mass of brass sample dissolved = 4.04 g

Absorbance of brass solution = 0.44

---

**OCR**  
Oxford Cambridge and RSA

**Copyright Information**

OCR is committed to seeking permission to reproduce all third-party content that it uses in its assessment materials. OCR has attempted to identify and contact all copyright holders whose work is used in this paper. To avoid the issue of disclosure of answer-related information to candidates, all copyright acknowledgements are reproduced in the OCR Copyright Acknowledgements Booklet. This is produced for each series of examinations and is freely available to download from our public website ([www.ocr.org.uk](http://www.ocr.org.uk)) after the live examination series.

If OCR has unwittingly failed to correctly acknowledge or clear any third-party content in this assessment material, OCR will be happy to correct its mistake at the earliest possible opportunity.

For queries or further information please contact The OCR Copyright Team, The Triangle Building, Shaftesbury Road, Cambridge CB2 8EA.

OCR is part of Cambridge University Press & Assessment, which is itself a department of the University of Cambridge.